

This Page Is Inserted by IFW Operations  
and is not a part of the Official Record

## BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images may include (but are not limited to):

- BLACK BORDERS
- TEXT CUT OFF AT TOP, BOTTOM OR SIDES
- FADED TEXT
- ILLEGIBLE TEXT
- SKEWED/SLANTED IMAGES
- COLORED PHOTOS
- BLACK OR VERY BLACK AND WHITE DARK PHOTOS
- GRAY SCALE DOCUMENTS

**IMAGES ARE BEST AVAILABLE COPY.**

**As rescanning documents *will not* correct images,  
please do not report the images to the  
Image Problem Mailbox.**




M-1f



The  
Patent  
Office

09/508133  
PCT B 99/01790

GB 99/01790

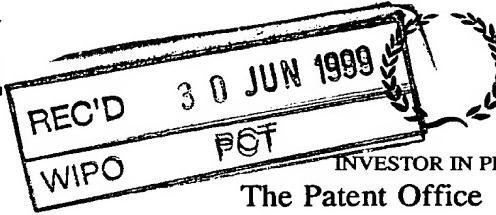

INVESTOR IN PEOPLE

The Patent Office

Concept House

Cardiff Road

Newport

South Wales

NP9 1RH

EJV

## PRIORITY DOCUMENT

SUBMITTED OR TRANSMITTED  
BUT NOT IN COMPLIANCE  
WITH RULE 17.1(a) OR (b)

I, the undersigned, being an officer duly authorised in accordance with Section 74(1) and (4) of the Deregulation & Contracting Out Act 1994, to sign and issue certificates on behalf of the Comptroller-General, hereby certify that annexed hereto is a true copy of the documents as originally filed in connection with the patent application identified therein.

I also certify that by virtue of an assignment registered under the Patents Act 1977, the application is now proceeding in the name as substituted.

In accordance with the Patents (Companies Re-registration) Rules 1982, if a company named in this certificate and any accompanying documents has re-registered under the Companies Act 1980 with the same name as that with which it was registered immediately before re-registration save for the substitution as, or inclusion as, the last part of the name of the words "public limited company" or their equivalents in Welsh, references to the name of the company in this certificate and any accompanying documents shall be treated as references to the name with which it is so re-registered.

In accordance with the rules, the words "public limited company" may be replaced by p.l.c., plc, P.L.C. or PLC.

Re-registration under the Companies Act does not constitute a new legal entity but merely subject the company to certain additional company law rules.

Signed

Dated 2 March 1999

THIS PAGE BLANK (USPTO)



The Patent  
Office

GB9819286.7

By virtue of a direction given under Section 30 of The Patents Act 1977, the application is proceeding in the name of

Michael Cole  
Poplar Farm  
Molesford  
WOODBRIDGE  
Suffolk  
IP13 0AL

**THIS PAGE BLANK (USPTO)**

25 SEP 1998

07SEP98 E387952-3 002819  
P01/7700 25.00 - 9819286.7

The Patent Office

**Request for grant of a patent**

(See the notes on the back of this form. You can also get an explanatory leaflet from the Patent Office to help you fill in this form)

THE PATENT OFFICE  
E  
05 SEP 1998

RECEIVED BY POST

Cardiff Road  
Newport  
Gwent NP9 1RH

1. Your reference

C319.00/G

2. Patent application number

(The Patent Office will fill in this part)

9819286.7

3. Full name, address and postcode of the or of each applicant (*underline all surnames*)SECTION 32  
Genevac Ltd  
The Sovereign Centre  
Farthing Road  
Ipswich IP1 5ARPatents ADP number (*if you know it*)

If the applicant is a corporate body, give the country/state of its incorporation

United Kingdom

07684743501

4. Title of the invention

Control of weight during evaporation of samples

5. Name of your agent (*if you have one*)Keith W Nash & Co  
90-92 Regent Street  
Cambridge  
CB2 1DP"Address for service" in the United Kingdom to which all correspondence should be sent (*including the postcode*)Patents ADP number (*if you know it*)

1206001 ✓

6. If you are declaring priority from one or more earlier patent applications, give the country and the date of filing of the or of each of these earlier applications and (*if you know it*) the or each application numberCountry      Priority application number  
(*if you know it*)      Date of filing  
(*day / month / year*)

7. If this application is divided or otherwise derived from an earlier UK application, give the number and the filing date of the earlier application

Number of earlier application      Date of filing  
(*day / month / year*)8. Is a statement of inventorship and of right to grant of a patent required in support of this request? (*Answer 'Yes' if:*

- a) any applicant named in part 3 is not an inventor, or
- b) there is an inventor who is not named as an applicant, or
- c) any named applicant is a corporate body.

yes c.)

See note (d))

9. Enter the number of sheets for any of the following items you are filing with this form.  
Do not count copies of the same document

Continuation sheets of this form

Description	5
Claim(s)	
Abstract	
Drawing(s)	2

10. If you are also filing any of the following, state how many against each item.

Priority documents

Translations of priority documents

Statement of inventorship and right to grant of a patent (*Patents Form 7/77*)

Request for preliminary examination and search (*Patents Form 9/77*)

Request for substantive examination  
(*Patents Form 10/77*)

Any other documents  
(please specify)

11.

I/We request the grant of a patent on the basis of this application.

Signature

Date

04/09/1998

KEITH W NASH & CO, AGENTS

12. Name and daytime telephone number of person to contact in the United Kingdom

C. Dieter Benziger (01223) 355477

#### Warning

After an application for a patent has been filed, the Comptroller of the Patent Office will consider whether publication or communication of the invention should be prohibited or restricted under Section 22 of the Patents Act 1977. You will be informed if it is necessary to prohibit or restrict your invention in this way. Furthermore, if you live in the United Kingdom, Section 23 of the Patents Act 1977 stops you from applying for a patent abroad without first getting written permission from the Patent Office unless an application has been filed at least 6 weeks beforehand in the United Kingdom for a patent for the same invention and either no direction prohibiting publication or communication has been given, or any such direction has been revoked.

#### Notes

- a) If you need help to fill in this form or you have any questions, please contact the Patent Office on 0645 500505.
- b) Write your answers in capital letters using black ink or you may type them.
- c) If there is not enough space for all the relevant details on any part of this form, please continue on a separate sheet of paper and write "see continuation sheet" in the relevant part(s). Any continuation sheet should be attached to this form.
- d) If you have answered 'Yes' Patents Form 7/77 will need to be filed.
- e) Once you have filled in the form you must remember to sign and date it.
- f) For details of the fee and ways to pay please contact the Patent Office.

C319.00/G

Title: Control of weight during evaporation of samples

 of Invention

This invention relates to a method of and apparatus for controlling the weight of samples dissolved or suspended in a liquid while they are evaporating in a vacuum. It is particularly applicable to samples in centrifugal evaporators.

Background to the Invention

Samples to be evaporated in centrifugal evaporators are usually held in glass or plastic tubes or, sometimes, in a large number of small wells in plastic blocks. The sample holders are mounted upon a rotating assembly and spun at relatively high speed so that a considerable centrifugal force is applied to them in an outward direction, which forces the liquid to the lower part of the sample tubes and prevents any frothing or spitting of the liquid out of the sample tubes when a vacuum is applied. The spinning samples are held in a vacuum-tight chamber (referred to herein as a "chamber") which is connected to a vacuum pumping device.

Centrifugal evaporators of this type are well known and many types are available commercially. One problem from which such evaporators suffer, is that it is very difficult to obtain a desired continuous read-out of the weight of the sample in the holders as the liquid is being evaporated, since the holders are being spun at a high speed, typically at about 1400 r.p.m. The possibility has been considered of continuously weighing the whole evaporator during spinning. However, this involves measuring a total weight of the order of 50 kg to an accuracy of about 1 gm, which is a very demanding task.

Another problem arises when evaporation needs to take place simultaneously for different

THIS PAGE BLANK (USPTO)

solvents, or solvent mixtures of differing compositions, in which the samples are dissolved or suspended. In this situation those samples which are held in the more volatile solvents or mixtures will evaporate faster than the ones held in the less volatile solvents, and this can lead to an excessive imbalance in the rotating assembly, and consequent unwanted vibrations. This would also mitigate against the possibility of weighing the whole evaporator.

In most centrifugal evaporator machines such unwanted vibrations are arranged to trip an out-of-balance sensor to thereby stop the machine, but in machines without a sensor the vibrations can cause damage to the machine and even to the samples. Sometimes the vibration problem can be overcome by careful loading of the evaporator, or by stopping the process from time to time and rebalancing the load by adding liquid to empty samples or by rearranging the samples in the rotating assembly. Both these methods are tedious and time consuming.

It is an object of the present invention to enable the weight of a sample in a centrifugal evaporator to be continuously and accurately measured during evaporation.

It is another object of the invention to enable the operation of a centrifugal evaporator to continue despite a considerable imbalance of forces.

#### Summary of the invention

According to one aspect of the invention there is provided a method of weighing a sample dissolved or suspended in a liquid during evaporation in a chamber of a centrifugal evaporator, comprising the steps of attaching a force measuring transducer between a sample holder and a rotor, obtaining force signals from the transducer while the rotor is spinning, and converting the signals into corresponding weights by taking account of a factor for the centrifugal force exerted on the sample holder.

**THIS PAGE BLANK (USPTO)**

The invention also extends to apparatus for measuring the weight of a sample as aforesaid, comprising a vacuum chamber, a sample holder connected to a rotor therein, a force measuring transducer attached between the holder and the rotor, and means for transmitting signals from the transducer to a reception means outside the chamber, and [REDACTED] for converting the signal to a value proportional to the weight of the sample holder.

The weight of the sample can then be continuously calculated by taking account of the centrifugal force and deducting the known weight of the holder.

The transducer is preferably a load cell; alternatively a strain gauge or position sensor may be utilised.

According to another aspect of the invention in the processing of samples in a centrifugal evaporator in which the samples are dissolved or suspended in liquids of differing volatility, any imbalance caused during spinning of the rotor and resulting in unwanted vibration is at least partially compensated for by associating with the rotor an automatic balancing unit.

Preferably the said unit comprises a bearing raceway incorporating a plurality of ball bearings which in use automatically distribute themselves to counteract the said imbalance.

#### Brief description of the drawings

The invention will now be described, by way of example only, with reference to the accompanying drawings in which:

Figure 1 is a schematic side view of a centrifugal evaporation system incorporating a force measuring transducer in accordance with the invention; and

THIS PAGE BLANK (USPTO)



Figure 2 is a perspective view of a disassembled automatic balancing unit associated with the rotor of Figure 1.

Detailed description

Referring first to Figure 1, this depicts a complete centrifugal evaporating system which is more fully described in copending Application No. 9803684.1 (ref: C311/G). The subject matter of said Application is therefore hereby entirely incorporated by reference into the present Application.

Sample holders in the form of blocks 4 are connected to pivots 13 such that the samples are held vertical during loading into the chamber 14. When the rotor shaft 5 is rotated the blocks swing out to a horizontal position due to centrifugal force.

Attached between each blocks and the shaft is a load cell 17 which records the horizontal force on the block which, when the rotor is spinning, will be proportional to the weight of the sample plus the sample holding assembly. Since the latter is constant the sample weight can readily be obtained. Of course, the apparent weight will be exaggerated by a factor due to the centrifugal force, but this factor will remain constant for a given rotor speed.

The power supply for the load cell may be provided in the same way as that for the direct temperature measurement described in the above mentioned Application No. 9803684.1. The weight signal is processed and transmitted to a receiver and decoder outside the chamber in a similar manner but via a separate channel.

In this way experiments have shown that one can measure weights of samples in holders weighing up to 1200 gm to an accuracy of better than 1 gm. This measurement allows accurate assessment of the correct moment at which to switch off sample heat, and of the point at which the samples are completely dry.

Figure 2 shows a proprietary automatic balancing unit which is preferably fitted to the

**THIS PAGE BLANK (USPTO)**

rotor shaft 5 as close as possible to the sample blocks 4. Here vibration caused by rotor imbalance is likely to occur, when solvents of different volatility are used for the samples.

The unit may be the Auto-Balancing unit produced by the bearing manufacturing company SKF.

As shown in Figure 2, the unit comprises inner and outer raceways 20 and 22 between which a number of loose ball bearings 24 are freely movable. The ball bearings distribute themselves automatically to counteract the imbalance in the rotor shaft 5.

**THIS PAGE BLANK (USPTO)**

# CENTRIFUGAL EVAPORATION SYSTEM

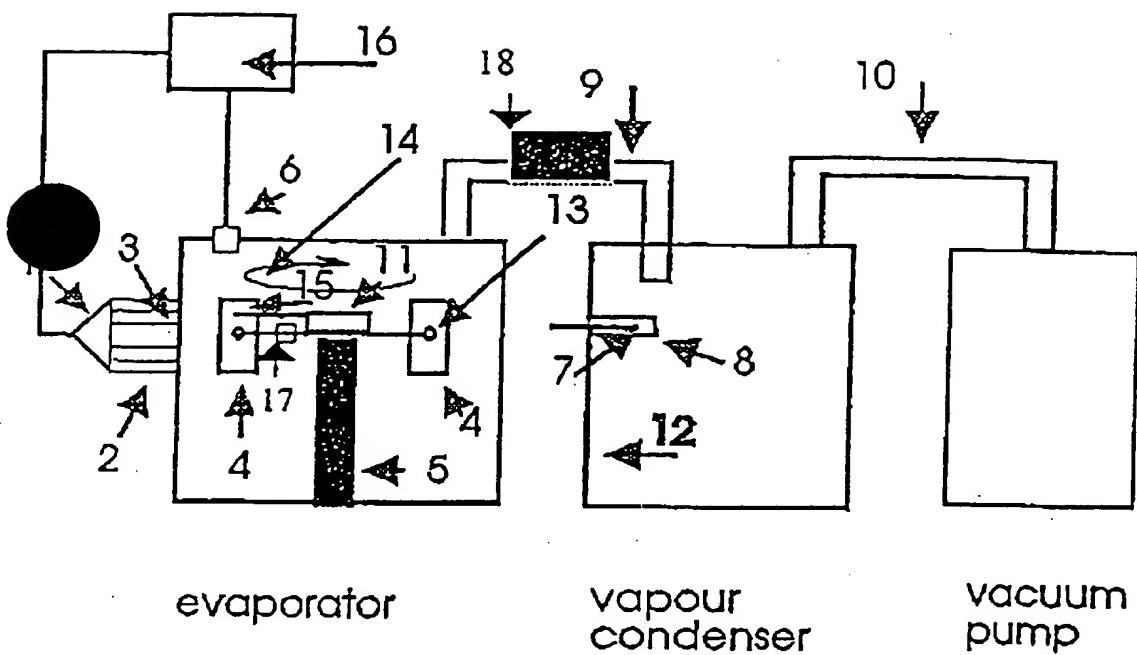


FIGURE 1

**THIS PAGE BLANK (USPTO)**

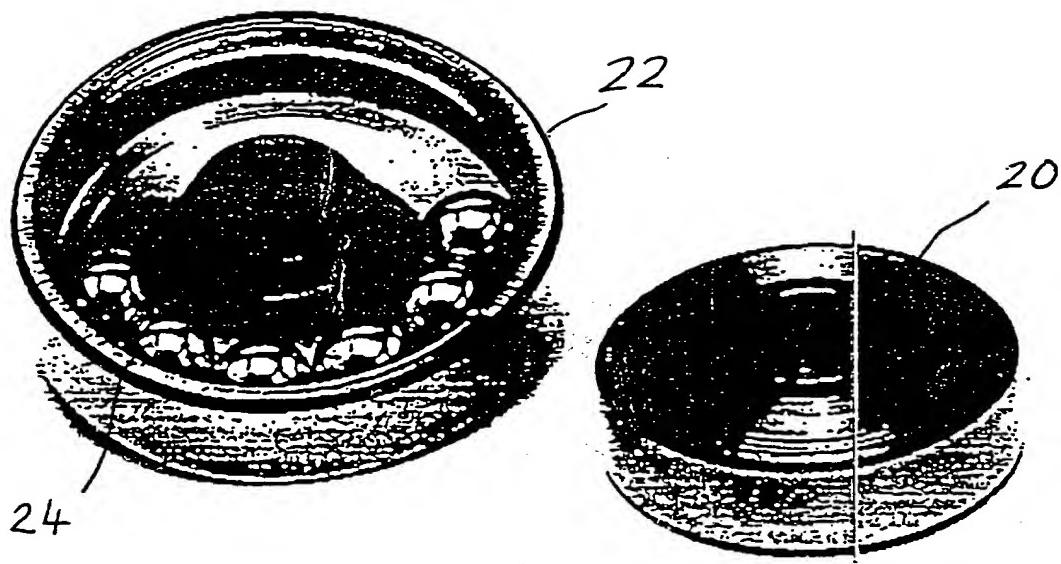


FIGURE 2

**THIS PAGE BLANK (USPTO)**